

Typenex® TypeSafe™ Segment Sampling Device

transfusion process steps:

Patient
Identification

Sample
Identification

Blood Bank
Component

Blood
Unit Issue

Recipient
Verification

Lab Supply
Basics

- TypeSafe allows for segment puncture without scissors for cross-matching procedures in the lab
- The large cannula allows cells to pass through with ease
- Top and bottom of device are easily identified by technicians in the blood bank
- Ribbed exterior allows for easy handling with gloved hands
- Fits into 10 and 12mL tube sizes
- Wide mouth fits most standard segment sizes

wide
mouth


ribbed
exterior

tapered
mid-section

large cannula



product is available in these colors:

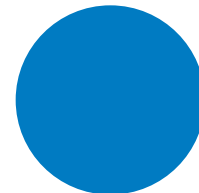

**4R4750
BLUE**

- **Highlights:** FDA cleared
- **Indications:** In-vitro diagnostic, single use
- **Packaging:** 1,000 units/case
(2 bags/case; 500 units/bag)

Enable safety in your laboratory with a Typenex lab supplies solution. Contact us for more information.

Phone | 866 897 3639
Fax | 866 897 3630
www.typenex.com